

Preoperative liver functional volumetry performed by 3D-99mTc-GSA scintigraphy/vascular fusion imaging using SYNAPSE VINCENT: a preliminary study

Hiroshi Yoshida¹, Hiroshi Makino¹, Tadashi Yokoyama¹, Hiroshi Maruyama¹, Atsushi Hirakata¹, Junji Ueda¹, Yasuhiro Mamada², Nobuhiko Taniai², Eiji Uchida²

¹Department of Surgery, Nippon Medical School Tama Nagayama Hospital, Tokyo 206-8512, Japan.

²Department of Surgery, Nippon Medical School, Tokyo 113-8603, Japan.

ABSTRACT

Aim: The present study was designed to evaluate the feasibility of preoperative liver functional volumetry performed by 3D-technetium-99m-diethylenetriaminepentaacetic acid-galactosyl-human serum albumin (99mTc-GSA) scintigraphy/vascular fusion imaging using SYNAPSE VINCENT and to examine the discrepancy between conventional and functional volumetry. **Methods:** The study group comprised 15 patients who underwent preoperative 3-dimensional (3D)-99mTc-GSA scintigraphy/vascular fusion imaging using SYNAPSE VINCENT software before hepatectomy between July 2014 and August 2015. The diagnosis was hepatocellular carcinoma ($n = 4$), metastatic liver tumor ($n = 10$), or intrahepatic cholangiocarcinoma ($n = 1$). Right hepatectomy was performed in 2 patients, left hepatectomy in 3 patients, right posterior sectionectomy in 3 patients, segmentectomy in 2 patients, and partial hepatectomy in 4 patients. 99mTc-GSA scintigraphy and computed tomography (CT) were performed to construct 3D-99mTc-GSA scintigraphy/vascular fused images. The conventional volume ratio of the planned resection region without tumor (% CT), and the functional volume ratio of the planned resection region without tumor (% GSA) were calculated. The discrepancy ratio was calculated as follows: discrepancy ratio = $100 - \% \text{GSA} / \% \text{CT} \times 100$ (%). **Results:** The % GSA ($17.9 \pm 16.7\%$) was significantly lower than the % CT ($21.5 \pm 17.6\%$) ($P < 0.036$). In all except 2 patients, the % GSA was lower than the % CT. The discrepancy ratio ranged from -4% to 75% (median, 20.7%). **Conclusion:** 3D-99mTc-GSA scintigraphy/vascular fused images constructed using SYNAPSE VINCENT were useful for noninvasively performing functional liver volumetry in patients scheduled to undergo various patterns of hepatectomy. In planned resection regions without tumor, the functional volume ratio was about 20% lower than the conventional volume ratio.

Key words: Functional volumetry; 99m-diethylenetriaminepentaacetic acid-galactosyl-human serum albumin; SYNAPSE VINCENT; fusion image; 3-dimensional; computed tomography

Address for correspondence:

Dr. Hiroshi Yoshida, Department of Surgery, Nippon Medical School Tama Nagayama Hospital, 1-7-1 Nagayama, Tama-city, Tokyo 206-8512, Japan. E-mail: hiroshiy@nms.ac.jp

Received: 25-03-2016, **Accepted:** 29-03-2016




Hiroshi Yoshida, M.D., is graduated from Nippon Medical School (Bsc) obtained MD and subsequently entered the Department of Surgery in 1986. From 1988 to 1992, He attended Postgraduate course in Department of Surgery, Nippon Medical School and obtained PhD. In 2003, he was appointed assistant professor and, in 2005, associated professor of the Department of Surgery, Nippon Medical School. In 2011, he was appointed Chief and Clinical Professor of the Department of Surgery, Nippon Medical School Tama Nagayama Hospital. In 2016, he was appointed Director of the Nippon Medical School Tama Nagayama Hospital.

This is an open access article distributed under the terms of the Creative Commons Attribution-NonCommercial-ShareAlike 3.0 License, which allows others to remix, tweak, and build upon the work non-commercially, as long as the author is credited and the new creations are licensed under the identical terms.

For reprints contact: service@oaepublish.com

How to cite this article: Yoshida H, Makino H, Yokoyama T, Maruyama H, Hirakata A, Ueda J, Mamada Y, Taniai N, Uchida E. Preoperative liver functional volumetry performed by 3D-99mTc-GSA scintigraphy/vascular fusion imaging using SYNAPSE VINCENT: a preliminary study. *Hepatoma Res* 2016;2:187-92.

Access this article online

Website: http://www.hrjournal.net/	Quick Response Code 
DOI: 10.20517/2394-5079.2016.06	

INTRODUCTION

In liver surgery, preoperative treatment planning is defined in the context of the anatomical structure and the functional reserve of the liver. In patients who have damaged liver function or who are scheduled to undergo extended hepatectomy, the postoperative residual liver volume with adequate preservation of blood supply and drainage vessels is very important for the prevention of liver failure.^[1]

Virtual reality simulation on computed tomography (CT), magnetic resonance imaging (MRI), or ultrasonography (US) plays an important role in examining the anatomical structure of the liver. Recently 3-dimensional (3D) imaging techniques, such as 3D CT, 3D MRI, and 3D US, have been developed. To date, a number of methods and software systems have been developed for 3D surgical planning before liver surgery.^[2-6] Remnant liver volume can thus be determined (volumetry) before hepatectomy.

Technetium-99m-diethylenetriaminepentaacetic acid-galactosyl-human serum albumin (^{99m}Tc-GSA) is an analog ligand of asialoglycoprotein. ^{99m}Tc-GSA binds to the asialoglycoprotein receptor (ASGP-R) located specifically on hepatocytes. The ASGP-R concentration is helpful in evaluating the extent and progression of liver disease, so the hepatic uptake of ^{99m}Tc-GSA reflects the number of functioning hepatocytes.^[7-11] Before hepatic resection, however, it is difficult to correctly estimate the functional hepatocyte mass of the remnant liver.

Thus, ^{99m}Tc-GSA scintigraphy combined with single-photo emission computerized tomography (SPECT) and CT fused imaging has been used to estimate the future remnant liver function before hepatic resection.^[12-16] However, the planned resection region had to be set manually using a 2-dimensional CT display. It was difficult to estimate the local remnant liver function in detail.

The volume analyzer software SYNAPSE VINCENT (Fujifilm Medical Co., Tokyo, Japan) is a high-speed 3D image analysis system. Using previously captured CT or MRI, high-definition 3D images of organs and vessels can be reconstructed quickly.^[4] It has become easy to grasp the complex anatomical relations between the portal triad, hepatic veins, and local tumor by volume rendering. With this software, the surgeon can simulate various patterns of planned hepatectomy.^[4-6]

In the present study, we performed preoperative liver functional volumetry by 3D-^{99m}Tc-GSA scintigraphy/vascular fused imaging using SYNAPSE VINCENT and examined the discrepancy between conventional and functional volumetry.

METHODS

This retrospective study was approved by our institutional

review board; informed consent was obtained from all patients before ^{99m}Tc-GSA scintigraphy and CT. We retrospectively reviewed patients who had undergone liver surgery between July 2014 and August 2015 in the Department of Surgery of our hospital. Twenty-five patients preoperatively underwent 3D-^{99m}Tc-GSA scintigraphy/vascular fused imaging using SYNAPSE VINCENT.

The following exclusion criteria were applied: (1) a history of hepatectomy or portal embolization; (2) hilar cholangiocarcinoma with unilateral biliary drainage; and (3) hepatectomy for benign disease. A total of 15 patients (10 men and 5 women; age, 60 to 81 years; mean, 72.7 years) who agreed to undergo preoperative 3D-^{99m}Tc-GSA scintigraphy/vascular imaging using SYNAPSE VINCENT were studied. The diagnosis was hepatocellular carcinoma in 4 patients, metastatic liver tumor in 10 patients, and intrahepatic cholangiocarcinoma in 1 patient. Based on the Brisbane 2000 classification criteria,^[17] right hepatectomy was performed in 2 patients, left hepatectomy in 3 patients, right posterior sectionectomy in 3 patients, segmentectomy in 2 patients, and partial hepatectomy in 5 patients. Background of the liver was liver cirrhosis in 2 patients, chronic hepatitis in 1 patient, and normal liver in 12 patients. The planned resection region of the liver could be similarly resected in all patients.

3D-vascular imaging using SYNAPSE VINCENT

Preoperative enhanced CT was performed with a 64-multidetector-row CT scanner (Discovery CT 750 HD, GE Healthcare Japan, Co. Ltd., Tokyo, Japan) at 0.625-mm intervals. Four-phase contrast-enhanced CT was performed 30, 60, 90, and 180 s after initiating the injection of contrast media to obtain hepatic arterial, portal venous, hepatic venous, and equilibrium phase images, respectively. A total of 100 mL of nonionic contrast material containing 370 mg of iodine per milliliter (Iopamidol, Bayer Yakuhin, Osaka, Japan) was injected intravenously at a rate of 3.3 mL/s using an automatic power injector. With the use of a workstation, a routine preoperative CT workup was performed in the axial and coronal imaging planes. The data were obtained in Digital Imaging and Communications in Medicine format and transmitted to a workstation running SYNAPSE VINCENT.

Liver Analysis Application uses Dynamic-CT imaging of the liver. After data for the hepatic arterial, portal venous, and hepatic venous phases are obtained, operative simulation by 3D-vascular images is performed. After selection of the portal venous branch of the planned resection area, surgical simulations can be displayed. This system also can calculate the total liver volume, tumor volume, and volume of planned resection region.^[4]

^{99m}Tc-GSA scintigraphy

A bolus of 1 mL of ^{99m}Tc-GSA (185 MBq, Nihon Medi-physics Co. Ltd., Nishinomiya, Japan) was intravenously injected

into an antecubital vein. Images were obtained with the patient in the supine position, using a gamma camera over a large field of view in which a high-resolution, all-purpose parallel-hole collimator (Infinia; GE Healthcare Japan Co. Ltd., Tokyo, Japan) was centered over the liver and precordium. Computer acquisition of gamma camera data was initiated simultaneously with injection of ^{99m}Tc -GSA and stopped at 30 min. Digital images (128×128 pixels) were acquired in the byte mode at a rate of 2 frames/min for 20 min after the injection. Hepatic SPECT data were obtained for 15 min after the end of the dynamic scintigraphic study.

3D- ^{99m}Tc -GSA scintigraphy/vascular fusion imaging using SYNAPSE VINCENT

Data obtained by ^{99m}Tc -GSA scintigraphy and CT imaging are composited by adjusting the axial and coronal images, and 3D- ^{99m}Tc -GSA scintigraphy/vascular fused images are constructed. The 3D-vascular images are used to select the portal venous branch to be resected and to calculate the extraction volume and ratio of the dominant region of the branch (conventional volumetry). In 3D- ^{99m}Tc -GSA scintigraphy, the functional volume and ratio of the same region are calculated (functional volumetry) concomitantly. Count-rates are displayed on the images, which can be saved.^[4]

Image analysis

In ^{99m}Tc -GSA scintigraphy, regions of interest (ROI) over the entire liver and heart are delineated. Time-activity curves are generated for the ROI.

An index of clearance of ^{99m}Tc -GSA from the blood is calculated as the quotient of the radioactivity of the heart ROI 15 min after injection (H15) divided by the radioactivity of the heart ROI 3 min after injection (H3), ($\text{HH15} = \text{H15}/\text{H3}$). Hepatic uptake ratio of ^{99m}Tc -GSA is calculated by dividing the radioactivity of the liver ROI at 15 min (L15) by the sum of H15 and L15 ($\text{LHL15} = \text{L15}/[\text{H15} + \text{L15}]$).^[10,11,18]

The 3D- ^{99m}Tc -GSA scintigraphy/vascular fusion images obtained using SYNAPSE VINCENT are used to calculate the total liver volume without tumor, the conventional volume ratio of the planned resection region without tumor (% CT), and the functional volume ratio of the planned resection region without tumor (% GSA). The discrepancy ratio is calculated as follows: discrepancy ratio = $100 - \% \text{GSA}/\% \text{CT} \times 100$ (%).

Case 9

A 71-year-old woman was admitted with a diagnosis of hepatocellular carcinoma concomitant with liver cirrhosis due to autoimmune hepatitis. CT revealed hepatocellular carcinoma (hypervascular tumor) in Segment 6 [Figure 1]. The % CT and % GSA were 3.0% and 2.2% in partial hepatectomy, 12.2% and 11.7% in segmentectomy, 27.5% and 28.4% in

right posterior sectionectomy, and 63.4% and 68.5% in right hepatectomy, respectively. Segmentectomy was performed to treat the hepatocellular carcinoma.

Statistical analysis

Values are expressed as means \pm standard deviation (range). The Student's *t* test was used to compare differences between two groups. A *P* value of < 0.05 was significant. All analyses were performed with SPSS version 17 (SPSS, Chicago, IL).

RESULTS


The clinical characteristics of the 15 patients are summarized in Table 1. The mean values of HH15 and LHL15 were 0.64 ± 0.10 and 0.90 ± 0.06 , respectively. The % GSA ($17.9\% \pm 16.7\%$) was significantly lower than the % CT ($21.5\% \pm 17.6\%$) ($P < 0.036$). In all except 2 patients, the % GSA was lower than the % CT. The discrepancy ratio ranged from -4% to 75% (median, 20.7%).

DISCUSSION

Our study demonstrated that 3D- ^{99m}Tc -GSA scintigraphy/vascular fusion imaging using SYNAPSE VINCENT is useful for noninvasive functional liver volumetry in patients scheduled to undergo various patterns of hepatectomy.

Postoperative liver failure remains a life-threatening complication after hepatectomy. Conventionally, traditional liver function tests and CT volumetry have been used to evaluate patients before hepatic surgery.

The use of ^{99m}Tc -GSA scintigraphy to evaluate liver function was initially introduced by investigators in Japan. ^{99m}Tc -GSA binds to the ASGP-R located specifically on hepatocytes. The function of ASGP-R remains normal even in regenerating hepatocytes, and the ASGP-R density per hepatocyte is constant. Therefore the total amounts of ASGP-R are lower in cirrhotic liver patients than in normal liver patients according to the hepatocyte theory.^[19-23] We previously reported that, with progression of hepatic functional degeneration, ASGP-R density per hepatic volume decreases, especially in the right lobe.^[9] The hepatic accumulation of ^{99m}Tc -GSA thus reflects the functional liver volume.^[24] The hepatic uptake image of ^{99m}Tc -GSA at 15 min or later reflects the ASGP-R population.^[25] An index of clearance of ^{99m}Tc -GSA (HH15) is calculated as the quotient of the radioactivity of the heart ROI 15 min after injection (H15) divided by the radioactivity of the heart ROI 3 min after injection (H3), ($\text{HH15} = \text{H15}/\text{H3}$). Hepatic uptake ratio of ^{99m}Tc -GSA (LHL15) is calculated by dividing the radioactivity of the liver ROI at 15 min (L15) by the sum of H15 and L15 ($\text{LHL15} = \text{L15}/[\text{H15} + \text{L15}]$).^[10,11,18] HH15 and LHL15 reflect the hepatic function. Various studies of ^{99m}Tc -GSA have examined hepatic function.^[9,10,26-28]

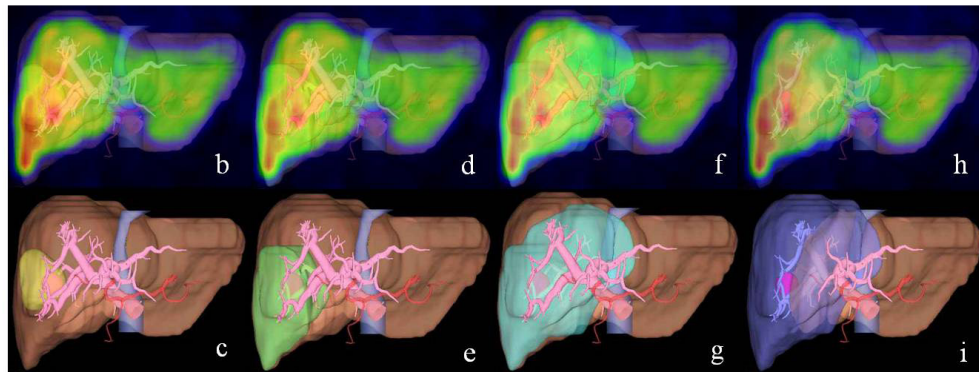

%GSA	2.2%	11.7%	28.4%	68.5%
%CT	3.0%	12.2%	27.5%	63.4%

Figure 1: Various patterns of % GSA and % CT on 3D-99mTc-GSA scintigraphy/vascular fusion images in case 11. (a) CT image (planned resection region); (b) partial hepatectomy; (c) partial hepatectomy; (d) segmentectomy (S6); (e) segmentectomy (S6); (f) right lateral sectionectomy; (g) right lateral sectionectomy; (h) right hepatectomy on 3D-99mTc-GSA scintigraphic image; (i) right hepatectomy on 3D-vascular image. % CT: conventional volume ratio of planned resection region without tumor; % GSA: functional volume ratio of planned resection region without tumor

Table 1: Discrepancy ratio between % CT and % GSA

Case	Age	Gender	Disease	Liver	Operation	HH15	LHL15	% CT	% GSA	Discrepancy ratio%
1	60	M	HCC	CH	Right hepatectomy	0.71	0.84	33.2	33.7	-2
2	76	M	Meta	Normal	Right hepatectomy	0.74	0.94	66.3	66.4	0
3	72	M	HCC	Normal	Left hepatectomy	0.44	0.94	42.7	19.0	56
4	65	F	ICC	Normal	Left hepatectomy	0.56	0.92	21.6	12.5	42
5	78	F	Meta	Normal	Left hepatectomy	0.62	0.90	32.9	29.1	12
6	68	M	Meta	Normal	Right posterior sectionectomy	0.57	0.96	28.9	25.5	12
7	74	M	Meta	Normal	Right posterior sectionectomy	0.70	0.86	27.5	23.5	15
8	77	M	Meta	Normal	Right posterior sectionectomy	0.70	0.92	25.5	23.2	9
9	71	F	HCC	LC	Segmentectomy (S6)	0.74	0.88	12.2	11.7	4
10	77	M	Meta	Normal	Segmentectomy (S8)	0.61	0.92	17.0	12.0	29
11	81	F	HCC	LC	Partial hepatectom (S8)	0.88	0.72	6.6	5.1	23
12	77	F	Meta	Normal	Partial hepatectomy (S6)	0.56	0.96	1.3	1.1	15
13	78	M	Meta	Normal	Partial hepatectomy (S6)	0.54	0.94	1.6	0.4	75
14	69	M	Meta	Normal	Partial hepatectomy (S8)	0.60	0.92	2.6	2.7	-4
15	68	M	Meta	Normal	Partial hepatectomy (S8)	0.68	0.88	2.9	2.2	24
AV	72.7					0.64	0.90	21.5	17.9	20.7
SD	5.6					0.10	0.06	17.6	16.7	21.4

HCC: hepatocellular carcinoma; Meta: metastatic liver tumor; ICC: intrahepatic cholangiocarcinoma; LC: liver cirrhosis; CH: chronic hepatitis; % CT: conventional volume ratio of planned resection region without tumor; % GSA: functional volume ratio of planned resection region without tumor; Discrepancy ratio: $(1 - \% \text{GSA} \% \text{CT}) \times 100$; AV: average value; SD: standard deviation

The SYNAPSE VINCENT system has already been tested in various surgical fields,^[29-32] and techniques for its use have been refined. The system has been greatly helpful in previsualizing intraoperative scenarios. In fact, the time required for intraoperative confirmatory US was considerably shortened after the introduction of this technology. This system is relatively easy to operate, allowing even a novice user to create 3D images.

3D-^{99m}Tc-GSA scintigraphy/vascular fusion images could also be easily created with the use of SYNAPSE VINCENT. It was possible to select the portal venous branch to be resected on the 3D-vascular images, and the conventional

volume could be compared with the functional volume of the planned resection region of the liver. The % GSA and % CT could also be calculated and compared.

In conclusion, 3D-^{99m}Tc-GSA scintigraphy/vascular fusion imaging performed with the use of SYNAPSE VINCENT is useful for non-invasive functional liver volumetry in patients scheduled to undergo various patterns of planned hepatectomy. There was a discrepancy between the results of conventional and functional volumetry. In the planned resection region without tumor, the functional volume ratio estimated with SYNAPSE VINCENT was about 20% lower than the conventional

